

## PRODUCT DATA SHEET

Standing: 2024-03-21

### LAUDA Proline PV 15 C

Clear-view thermostat 230 V; 50/60 Hz

Part Number: L001535

#### Features

- Thermostatic bath/circulator with latest microprocessor technology and many functions for programming and information and installed cooling coil
- Back-lit graphic LCD display with high resolution and different display modes
- Additional green LED display for temperature
- Input either via cursor keys, numeric soft keys or both. Additional Tmax key for overtemperature
- Command console can be detached and used as remote control
- EasyUse system for simple operation of the whole unit
- SelfCheck Assistant for system diagnosis
- Fully electronic continuous controller with PID action for internal & external control
- PowerAdapt system for the use of the maximum possible amount of heat permitted by the power supply system
- Low-level and adjustable over-temperature protection with acoustic alarm for use with flammable and non-flammable liquids
- LAUDA Varioflex pump (pressure) with 8 selectable levels
- Optically decoupled RS 232/485 interface integrated as a standard
- Option for upgrading with up to 2 interfaces (RS 232/485, Profibus, analogue or contact modules, Ethernet-USB module)
- Programmer with 150 temperature/time segments that can be separated into 5 programs
- Timer function for switching on the thermostat, entering the standby mode or running of programs
- Measuring chamber with big windows front and back
- Pump connectors on the side and in the back, installed bypass
- Bath vessel made of stainless steel (insulated, with handles and drain tap)



Working temperature min.  
30 °C



Working temperature max.  
230 °C

LAUDA DR. R. WOBSEY GMBH & CO. KG  
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0  
info@lauda.de • www.lauda.de  
WEEE-Reg.-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen  
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:  
LAUDA DR. R. WOBSEY Verwaltungs-GmbH  
Sitz Lauda-Königshofen  
Registergericht Mannheim • HRB 560226

Geschäftsführer:  
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,  
Dr. Ralf Hermann, Dr. Marc Stricker  
Beirat: Dr. Gerhard Wobser

# PRODUCT DATA SHEET

Standing: 2024-03-21

LAUDA Proline PV 15 C

Clear-view thermostat 230 V; 50/60 Hz

Part Number: L001535

## Technical Features (according to DIN 12876)

Working temperature range	30 ... 230 °C
Working temperature range with external cooling	20 ... 230 °C
Operating temperature range	0 ... 230 °C
Ambient temperature range	5 ... 40 °C
Temperature stability	0.01 ± K
Heater power max.	3.6 kW
Current max.	16 A
Power consumption max.	3.7 kW
Pump Pressure max.	0,8 bar
Pump flow rate max. (pressure)	25 L/min
Bath volume min. / max.	11.0 / 15.0 L
Size of bath (W x D x H)	230 x 135 x 320 mm
Overall dimensions (WxDxH)	506 x 282 x 646 mm
Weight	31 kg
Power supply	230 V; 50/60 Hz
Power plug	Power cord with angled plug (CEE7/7)

Reserve technical changes

## Standard accessories

- 4 screw caps, 4 closing plugs
- 2 nipples 13 mm for pump connectors
- 2 nipples for cooling coil

LAUDA DR. R. WOBSEY GMBH & CO. KG  
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0  
info@lauda.de • www.lauda.de  
WEEE-Reg.-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen  
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:  
LAUDA DR. R. WOBSEY Verwaltungs-GmbH  
Sitz Lauda-Königshofen  
Registergericht Mannheim • HRB 560226

Geschäftsführer:  
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,  
Dr. Ralf Hermann, Dr. Marc Stricker  
Beirat: Dr. Gerhard Wobser


## PRODUCT DATA SHEET

Standing: 2024-03-21

LAUDA Proline PV 15 C

Clear-view thermostat 230 V; 50/60 Hz

Part Number: L001535

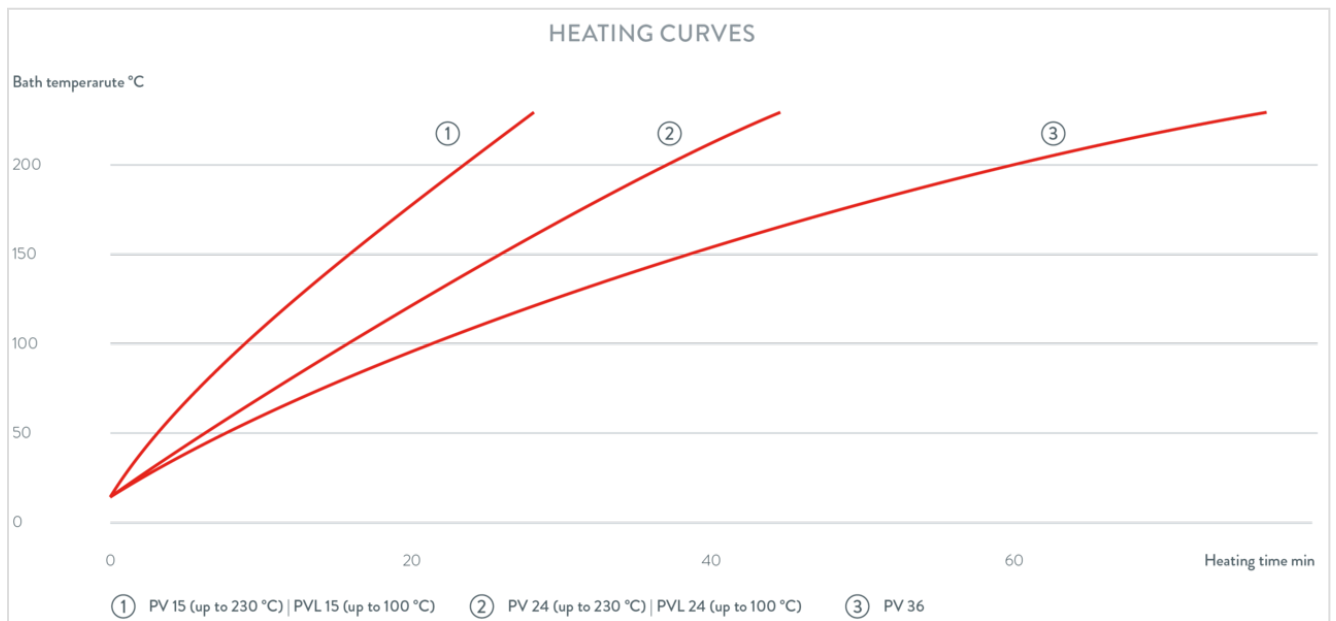

Reserve technical changes

LAUDA DR. R. WOBSEY GMBH & CO. KG  
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0  
info@lauda.de • www.lauda.de  
WEEE-Reg-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen  
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:  
LAUDA DR. R. WOBSEY Verwaltungs-GmbH  
Sitz Lauda-Königshofen  
Registergericht Mannheim • HRB 560226

Geschäftsführer:  
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,  
Dr. Ralf Hermann, Dr. Marc Stricker  
Beirat: Dr. Gerhard Wobser


## PRODUCT DATA SHEET

Standing: 2024-03-21

LAUDA Proline PV 15 C

Clear-view thermostat 230 V; 50/60 Hz

Part Number: L001535

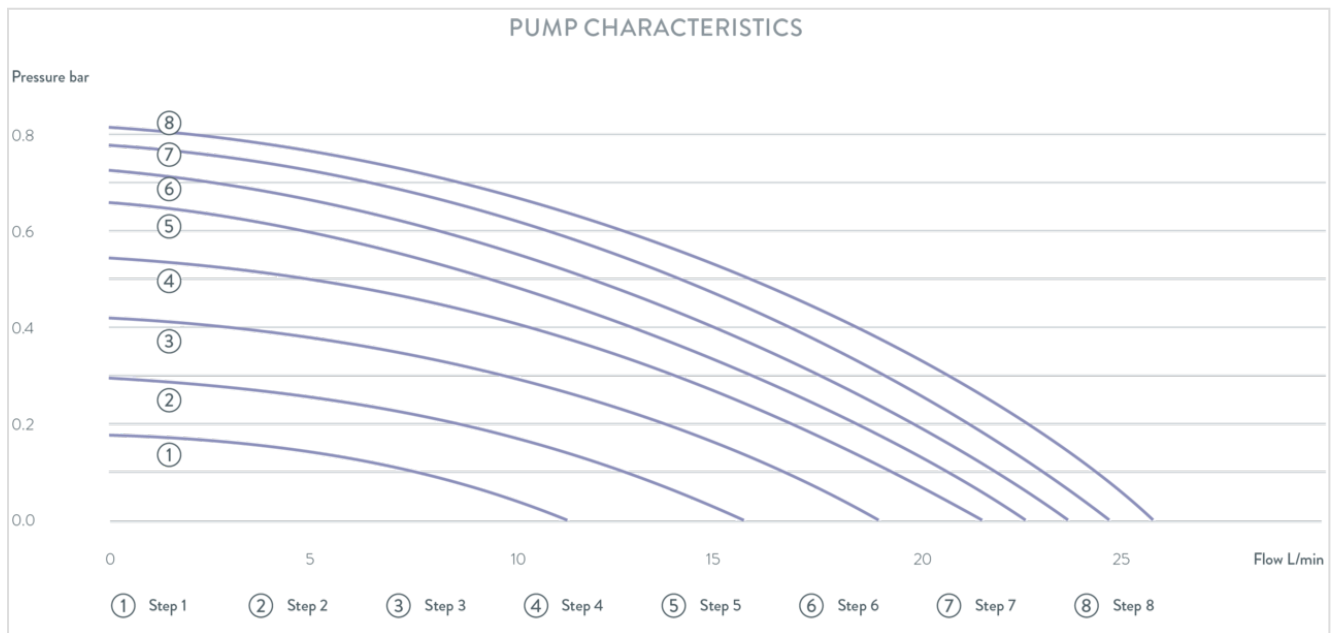

Reserve technical changes

LAUDA DR. R. WOBSEY GMBH & CO. KG  
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0  
info@lauda.de • www.lauda.de  
WEEE-Reg-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen  
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:  
LAUDA DR. R. WOBSEY Verwaltungs-GmbH  
Sitz Lauda-Königshofen  
Registergericht Mannheim • HRB 560226

Geschäftsführer:  
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,  
Dr. Ralf Hermann, Dr. Marc Stricker  
Beirat: Dr. Gerhard Wobser