

PRODUCT DATA SHEET

Standing: 2024-01-22

LAUDA ECO Silver

Immersion thermostat 230 V; 50/60 Hz

Part Number: L001076

Features

- Immersion thermostat with latest microprocessor technology
- Monochrome graphic LCD display for simultaneous indication of actual and set value
- User friendly menu navigation in plain language
- Easy input via cursor and soft keys. Additional Tmax key for overtemperature protection
- Fully electronic continuous controller with PID action
- Safety class III for operation with flammable and non-flammable liquids. Over-temperature cut-out adjustable via menu
- Vario pump with six adjustable performance levels
- Easy control of the flow rate between internal and external circulation during operation without contact with the bath
- USB interface as standard
- Upgradeable with an interface module (analogue module, contact module, RS 232/485 module, Profibus, Ethernet-USB module)
- Upgradeable with Pt 100/LiBus module for external control and remote control via Command console
- Integrated programmer, 1 programme with max. 20 segments
- With screw clamp for bath vessels with a wall thickness of up to 25 mm and an immersion depth of min. 150 mm



Reserve technical changes



Working temperature min.
20 °C



Working temperature max.
200 °C

LAUDA DR. R. WOBSEY GMBH & CO. KG
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0
info@lauda.de • www.lauda.de
WEEE-Reg-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:
LAUDA DR. R. WOBSEY Verwaltungs-GmbH
Sitz Lauda-Königshofen
Registergericht Mannheim • HRB 560226

Geschäftsführer:
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,
Dr. Ralf Hermann, Dr. Marc Stricker
Beirat: Dr. Gerhard Wobser

PRODUCT DATA SHEET

Standing: 2024-01-22

LAUDA ECO Silver

Immersion thermostat 230 V; 50/60 Hz

Part Number: L001076

Technical Features (according to DIN 12876)

Working temperature range	20 ... 200 °C
Ambient temperature range	5 ... 40 °C
Temperature stability	0.01 ± K
Heater power max.	2 kW
Power consumption max.	2.1 kW
Current max.	10 A
Pump Pressure max.	0,6 bar
Pump flow rate max. (pressure)	22 L/min
Overall dimensions (WxDxH)	130 x 135 x 325 mm
Weight	3 kg
Power supply	230 V; 50/60 Hz
Power plug	Power cord with angled plug (CEE7/7)

Standard accessories

- 1 screw clamp

LAUDA DR. R. WOBSEY GMBH & CO. KG
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0
info@lauda.de • www.lauda.de
WEEE-Reg.-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:
LAUDA DR. R. WOBSEY Verwaltungs-GmbH
Sitz Lauda-Königshofen
Registergericht Mannheim • HRB 560226

Geschäftsführer:
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,
Dr. Ralf Hermann, Dr. Marc Stricker
Beirat: Dr. Gerhard Wobser


PRODUCT DATA SHEET

Standing: 2024-01-22

LAUDA ECO Silver

Immersion thermostat 230 V; 50/60 Hz

Part Number: L001076

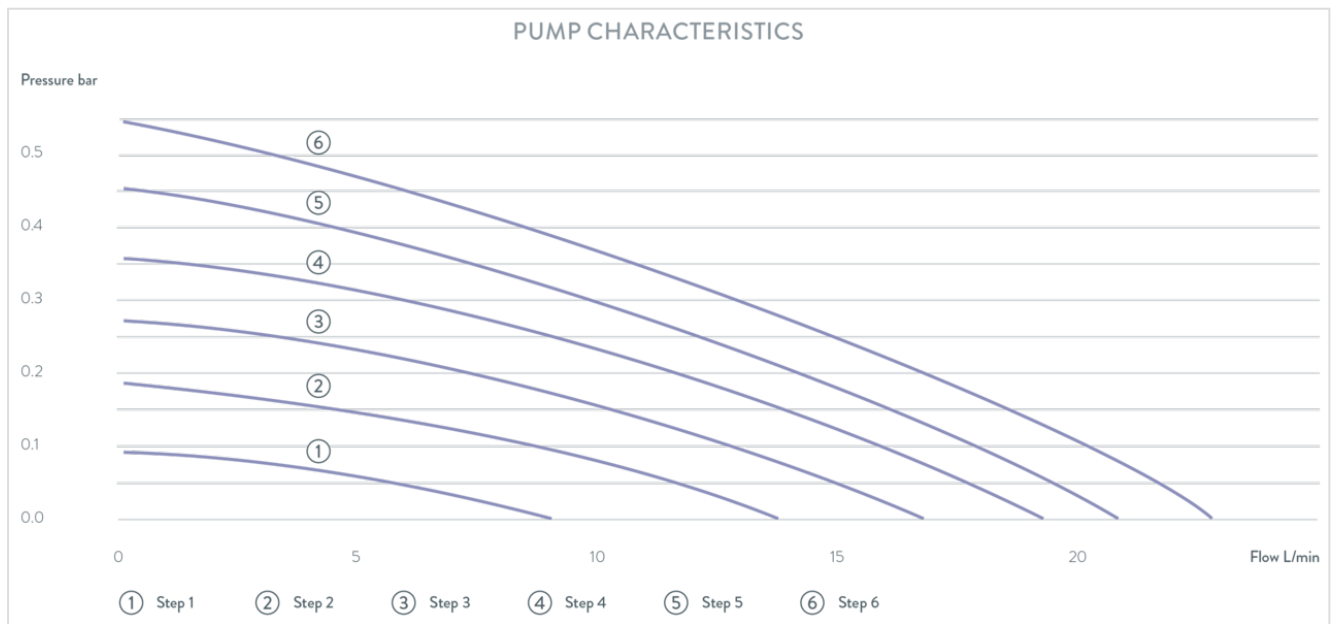

Reserve technical changes

LAUDA DR. R. WOBSEY GMBH & CO. KG
Laudaplatz 1 • 97922 Lauda-Königshofen • DE

T + 49 (0) 9343 503-0
info@lauda.de • www.lauda.de
WEEE-Reg-Nr.: DE 66 42 40 57

Kommanditgesellschaft: Sitz Lauda-Königshofen
Registergericht Mannheim • HRA 560069

Persönlich haftende Gesellschafterin:
LAUDA DR. R. WOBSEY Verwaltungs-GmbH
Sitz Lauda-Königshofen
Registergericht Mannheim • HRB 560226

Geschäftsführer:
Dr. Gunther Wobser (Vors.), Dr. Mario Englert,
Dr. Ralf Hermann, Dr. Marc Stricker
Beirat: Dr. Gerhard Wobser